

> Back To Main Page  
Paragon Medical Supply



**COVIDIEN**

*positive results for life™*

WarmTouch™ Convective Warming System.  
Maintaining normothermia has never been so easy.



## Covidien offers the **right solution**

The WarmTouch™ convective warming system provides veterinary clinicians with an easy solution to help maintain normothermia. The WarmTouch™ convective warming unit has an intuitive screen showing the temperature and status of each of the 5 temperature settings. The soft but strong WarmTouch™ blanket range is tailored for most surgical procedures, offering tear resistance, even airflow and patient comfort.



# The innovative WarmTouch™ Convective Warming Unit

## **EASY TO USE**

New intuitive LCD screen and menu function provide all information at a glance.

## **EASY TO CONTROL**

5 different settings, including ambient and boost, for more flexibility and precision.

## **EASY TO HANDLE**

The low weight makes it easy to transport and attach to a cart, bed rail or an IV pole.

## **EASY TO SERVICE**

The onscreen filter expiration warning makes it easy to maintain and service.




# WarmTouch™ Convective Warming Blankets: A distinctive design to optimize warming efficacy



The soft but strong WarmTouch quilted blankets meet a wide range of needs to offer the optimal solution in easy and effective patient warming.

Made of a strong two-ply material, they resist tears, punctures and fluids, yet are soft and flexible to allow custom draping and comfort for your patients.

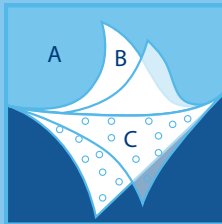

- A** The fabric outer layer eliminates hot plastic surfaces that could make contact with the patient.
- B** Layered construction adds strength, so even minor punctures will not expand and cause the blanket to deflate.
- C** Specially positioned holes direct the airflow evenly onto the patient.

# Warming with precision, handling with flexibility

A full range of blankets to accommodate a range of patients and procedures.



**WARMTOUCH™  
FULL BODY  
BLANKET**

---

Easy access to chest tubes, arteries and pulse checking.  
Also available for pediatric.



**WARMTOUCH™  
UPPER BODY  
BLANKET**

---

Dual inlets for added flexibility.



**WARMTOUCH™  
LOWER BODY  
BLANKET**

---

Precise, localized warming with adjustable tape positioning.



**WARMTOUCH™  
CARDIAC BLANKET  
STERILE**

---

Unique design to minimize air distribution at patient's feet.



**WARMTOUCH™  
PEDIATRIC  
SURGICAL BLANKET**

---

For over and under body use. Clear plastic head cover to maintain head temperature.

# WarmTouch™ Convective Warming System

## ordering information and specifications

| DESCRIPTION                  | SIZE  | QUANTITY | CATALOG NUMBER |
|------------------------------|---|----------|----------------|
| WarmTouch™ 6000 Warming Unit | 24" x 17" x 12" (60 cm x 43 cm x 30 cm)<br>(hose fully collapsed, nozzle in storage position on unit) | 1/box    | 5016000        |

### WARMING UNIT

**Weight with Power Cord** 5.2kg (11.5 pounds)

### TRANSPORT CART (OPTIONAL)

**Weight** 3.1kg (6.8 pounds)  
**Height** 67.1cm (26.4 inches)  
**Width** 32.3cm (12.7 inches)  
**Depth** 38.6cm (15.2 inches)

### ELECTRICAL REQUIREMENTS ENVIRONMENTAL

**Power Requirements** 100 to 240 V AC  
 Max. current at 100 V = 8 A  
 Max. current at 240 V = 5 A

**Input Frequency** 50/60 Hz ±1 Hz

**Fuse (x2)** Littlefuse (mfr.) 0218010; 250 V, 10 A

**Temperature** 18°C to 28°C (64.4°F to 82.4°F)

### PERFORMANCE SPECIFICATIONS

**Maximum Contact Surface Temperature** 44.1°C (111.4°F)

**Average Time for Contact Surface Temperature to Rise from 23°C ±2°C to 37°C** 6 minutes

**Average Time for Temperature of Air Exiting the Hose to Rise from 23°C ±2°C to 37°C** < 1 minute

**Accuracy of Displayed Temperature** ±1.0°C (air entering hose)

**Automatic Temperature Stepdown (Boost to High Temperature)** After 45 minutes of continuous use, blower will step down from the Boost to High setting.

**Thermal Protection Threshold** Thermostat (internal): 49°C to 55°C (120°F to 131°F)

**Average Alarm Level** 56 dB @ 1 meter



## WarmTouch™ Convective Warming Blankets size and ordering information

| DESCRIPTION                                 | SIZE  | QUANTITY | CATALOG NUMBER |
|---|---|----------|----------------|
| WarmTouch™ Full-Body / Multi Access Blanket | 40" x 79" uninflated (102 x 165 cm)<br>33" x 74" inflated (84 x 188 cm) | 12/case  | 5030810        |
| WarmTouch™ Pediatric Blanket                | 40" x 79" uninflated (102 x 165 cm)<br>33" x 74" inflated (84 x 188 cm) | 12/case  | 5030840        |
| WarmTouch™ Pediatric Surgical Blanket       | 25" x 41" uninflated (63 x 104 cm)<br>22" x 35" inflated (56 x 89 cm)   | 12/case  | 5030850        |
| WarmTouch™ Cardiac Blanket                  | 40" x 65" uninflated (102 x 165 cm)<br>34" x 53" inflated (86 x 135 cm) | 12/case  | 5030860        |
| WarmTouch™ Upper Body Blanket               | 28" x 82" uninflated (71 x 208 cm)<br>24" x 72" inflated (61 x 183 cm)  | 12/case  | 5030870        |
| WarmTouch™ Lower Body Blanket               | 41" x 56" uninflated (104 x 142 cm)<br>34" x 52" inflated (86 x 132 cm) | 12/case  | 5030880        |

**WarmTouch™ Convective Warming System.**  
Maintaining normothermia has never been so easy.



**WarmTouch™**  
Lower Body Blanket



**WarmTouch™**  
Convective Warming Unit



**COVIDIEN**

6135 GUNBARREL AVENUE  
BOULDER, CO 80301  
800-635-5267

[WWW.COVIDIEN.COM](http://WWW.COVIDIEN.COM)

**IMPORTANT:** Please refer to the package insert for complete instructions, contraindications, warnings and precautions.

COVIDIEN, COVIDIEN with logo, Covidien logo and *positive results for life* are U.S. and internationally registered trademarks of Covidien AG. Other brands are trademarks of a Covidien company. © 2013 Covidien. 13-PM-0157 TM90047